

EXTRACTION FORCEPS





Extraction forceps from American Eagle Instruments® with weight-reduced, ergonomic handles and a smooth, matt surface are particularly easy to clean.


The anatomically shaped beaks of the standard extraction forceps are designed with an adhesive profile for a **secure grip** - they can be used with little effort. A Teflon disc in the instrument connection maintains ease of movement and prevents corrosion.

There are also some special pliers in the range:

Deep Grip: Deep-reaching extraction forceps with slightly sharpened forcep beaks that are elongated and tapered towards the tip. They allow subgingival placement and distribute the pressure over a larger area of the often already weakened crown and root structure.

Root and residual root extraction forceps: S-curved Witzel forceps can be used almost universally due to its shape. **Only one root forceps for all quadrants!** The extended beak of the pliers, which is tapered towards the tip, has a diamond coating to ensure a secure grip even with very fine root residues.

There are also three other variants available with extra sharpened jaws or blades that ensure a secure hold: for particularly demanding root residue removal deep in the alveolus (AESRFLO, AESRFUN, AESRFUP).

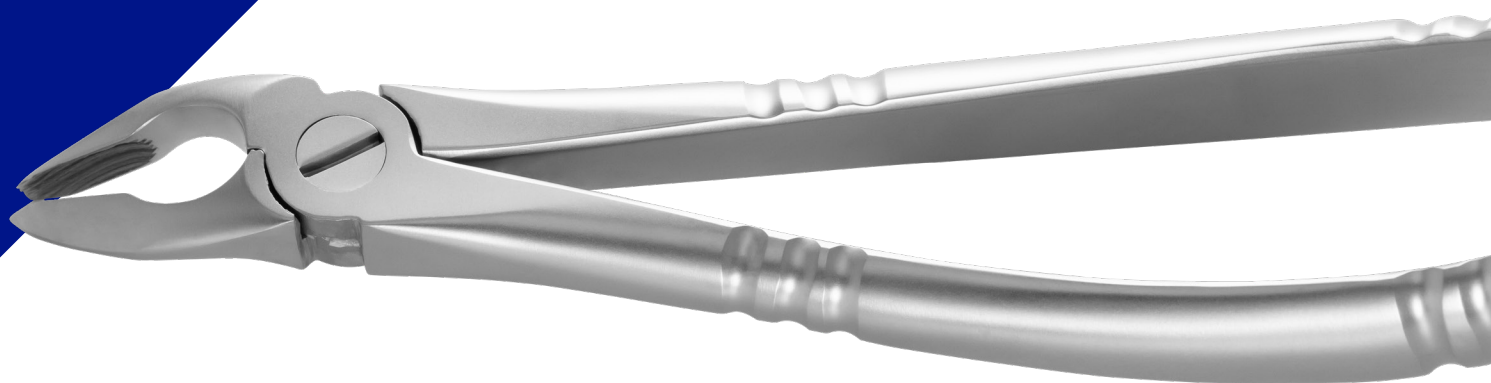

| | | | |
|----------|---|-----------|--|
| AESEFSET | Standard Extraction Forceps Set | | |
| | Standard Extraction Forceps Set, including 10 forceps: for upper incisors and canines, for upper premolars, for upper right molars, for upper left molars, for upper wisdom teeth, for upper roots, for lower canines and premolars, for lower molars both sides, for lower roots ad for lower wisdom teeth (AESEFCUI, AESEFCUP, AESEFCLI, AESEFCURM, AESEFCULM, AESEFCLM, AESEFCLR, AESEFCUR, AESEFCU3M, AESEFCL3M) | | |
| AESEFUI | Deep Grip Upper Incisor Extraction Forceps | AESEFCURM | Standard Extraction forceps #17 upper right molars |
| AESEFUC | Deep Grip Upper Canines and Premolars Extraction Forceps | AESEFCULM | Standard Extraction forceps #18 upper left molars |
| AESEFLCP | Deep Grip Lower Front, Canines and Premolars Extraction Forceps | AESEFCLM | Standard Extraction forceps #22 lower molars |
| AESEFR | Deep Grip Low Roots and Narrow Incisors Extraction Forceps | AESEFCLR | Standard Extraction forceps #33 lower roots |
| AESEFUM | Deep Grip Upper Molars and Wisdom Teeth Extraction Forceps | AESEFCUR | Standard Extraction forceps #51A upper roots |
| AESEFLM | Deep Grip Lower Molars and Wisdom Teeth Extraction Forceps | AESEFCU3M | Extracting Forceps #67A with Non-Slip Profile for Upper Wisdom Teeth on Both Sides |
| AESEFRW | Witzel Universal Root Extraction Forceps, Diamond Coated | AESEFCL3M | Extracting Forceps #79 with Non-Slip Profile for Lower Wisdom Teeth on Both Sides |
| AESEFCUI | Standard Extraction forceps #2 upper incisors and canines | AESRFLO | Root Forceps - lower |
| AESEFCUP | Standard Extraction forceps #7 upper premolars | AESRFUN | Root Forceps - universal |
| AESEFCLI | Standard Extraction forceps #13 lower incisors, Canines and premolars | AESRFUP | Root Forceps- upper |